

A photograph of two women in a retail or medical setting. The woman on the left has long brown hair and is wearing a white shirt with a measuring tape around her neck. The woman on the right has short dark hair and is wearing a grey top. They are both looking down at something out of frame with interest. In the background, there are shelves with various items, possibly medical supplies or retail products.




medi

Measuring chart set

Round-knitted
Flat-knitted



medi. I feel better.

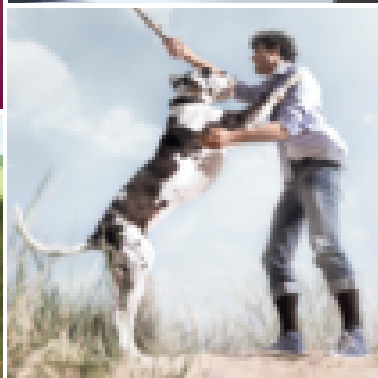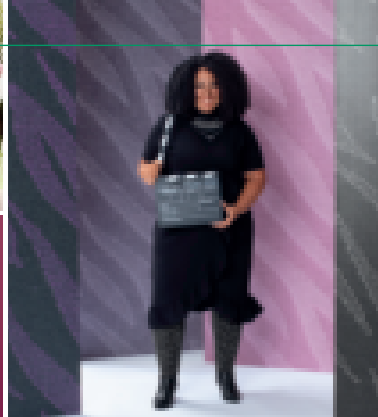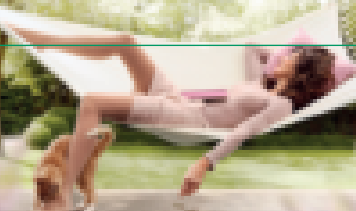



Driving the future in compression therapy

**A strong brand. A wide-ranging product range.
Convincing concepts.**

medi World of Compression: signpost to the future with more than 65 years' competence in the field of compression. Your partner for the highest demands on quality, innovation and design. Let's shape the future together and create positive brand experiences – for you and your clients!



medi

medi produces medical aids for a variety of conditions and life situations. But we don't see ourselves solely as a manufacturer. On the contrary, we handle our own research and development, technology, production, marketing and sales in-house – always keeping our focus on the market, its demands and variables, and the medi vision.

Headquarters
medi GmbH & Co. KG
Medicusstraße 1
95448 Bayreuth
Germany
T +49 921 912-0
F +49 921 912-781
export@medi.de
www.medi.de/en

Order hotline

T +49 921 912-0

Orders, delivery information, product information
and advice

Product service

T +49 921 912-0

F +49 921 912-781

Complaints processing standard, product suggestions

mediven® product features



Clima Comfort

Clima Comfort*

- Wicks moisture away from the skin's surface quickly
- Guarantees a high level of breathability



Clima Fresh

Integrated Clima Fresh function

- Antibacterial effect
- Thus inhibits odour-forming bacteria



Perfect Fit

Perfect Fit

- A very precise fit
- Medically effective pressure gradient
- Optimal care conditions in standard and made-to-measure version



Soft Elastic

Soft Elastic

- Noticeably more toe freedom with closed toe
- Increases wearing comfort
- Less pressure is exerted on the end of the foot

* Scientifically tested by French institute IFTH. The test results show that Clima Comfort technology from medi assures climate control between the skin and the stocking. With 11/12 points, medi achieved an outstanding overall score.

mediven[®] product features



Variety

A wide range of variants and special models



Anatomic Fit

Anatomical oblique borders as standard (distal and proximal)

- Unimpeded lymphatic drainage without circular constriction
- Optimum fit and better hold



Softness

Softness**

- Soft material for a pleasant wearing comfort
- Slides smoothly and gently over the skin



Antibacterial

Antibacterial**

- Knitted fabric with antibacterial effect
- Inhibits growth and nidation of malodorous bacteria



Easy donning

Easy donning***

- Liner and outer stocking with both 20 mmHg for easy donning and doffing
- Soft liner stocking eases donning of the outer stocking

** These features only apply to mediven ulcer liner stocking

*** These features only apply to mediven ulcer kit

mediven® product features

Hightech fibres for natural elasticity

- High double-stretch elasticity
- Optimal user comfort
- Easy to put on and take off
- Perfect, wrinkle-free fit
- Highly attractive transparency and refined feeling



cf. Ben Ammar, F. (Dermscan): Clinical evaluation of the skin compatibility of a medical device - application testing realized under dermatological control, Tunis 2012/2013.

Colours for greater quality of life

- Large selection of standard colours and the latest trend colours available.








All mediven compression stockings are latex-free.



Two steps to the perfect, individualised mediven® compression stocking

Step 1: Indication overview

			mediven elegance (1, 2)	mediven for men (1, 2)	mediven comfort (1, 2)	mediven plus (1, 2, 3)	mediven active (1, 2)	mediven forte (2, 3)	mediven harmony (2, 3)	mediven ulcer kit (2 x 20 mmHg)
Indication		Product (CCL)								
Chronic venous diseases according to CEAP										
C0s		No visible signs of a venous disease	●	●	●	●	●	●		
C1		Spider veins, telangiectasias, reticular varicose veins	●	●	●	●	●	●		
C2		Varicose veins with no signs of a CVI	●	●	●	●	●	●		
C3		Prevention of and therapy for venous oedema		●	●	●	●	●		
C4	 a)  b)	Prevention and therapy of skin changes C9a Pigmentation, eczema C9b Lipodermatosclerosis, atrophy blanche			●	●	●	●		
C5		Healed venous ulcers (UCV), Prevention of UCV recurrence				●		●		●
C6		Florid UCV								●

		Light material			Medium material		Strong material		
	Product (CCL)	mediven elegance (1, 2)	mediven for men (1, 2)	mediven comfort (1, 2)	mediven plus (1, 2, 3)	mediven active (1, 2)	mediven forte (2, 3)	mediven harmony (2, 3)	mediven ulcer kit (2 x 20 mmHg)
Indication									
Other									
Therapy for mixed venous ulcer									●
Other chronic venous diseases									
Initial phase following varicose therapy		●		●	●		●		
Functional venous insuffi ciency due to obesity							●		
Other oedemas									
Congestion complaints and oedema in pregnancy, postoperative reperfusion oedema, post-traumatic, postoperative, occupational oedema, medicinally induced oedema, oedema due to immobility			●	●	●	●	●		
Thromboembolism									
Superficial venous thrombosis		●	●	●	●	●	●		
Deep vein thrombosis					●		●		
Condition after thrombosis, post-thrombotic syndrome				●	●	●	●		
Thrombosis prophylaxis in mobile patients		●	●	●	●	●	●		
Arm vein thrombosis								●	
Other indications									
Inflammatory dermatoses of the legs		●	●	●	●	●	●		
Nausea, dizziness in pregnancy		●		●	●		●		

medi Tipp: Safe care for concomitant diseases.

Customers with chronic venous disease and simultaneously


→ mild to moderate peripheral arterial disease (PAD) and / or

→ diabetes mellitus / diabetic polyneuropathy (PNP)

have special requirements for a medical compression stocking.

mediven angio was particularly developed for this purpose – and its safety was proven in a clinical trial. Actively ask your customers about these concomitant diseases.

Step 2: Determine the indication and evaluate comorbidities

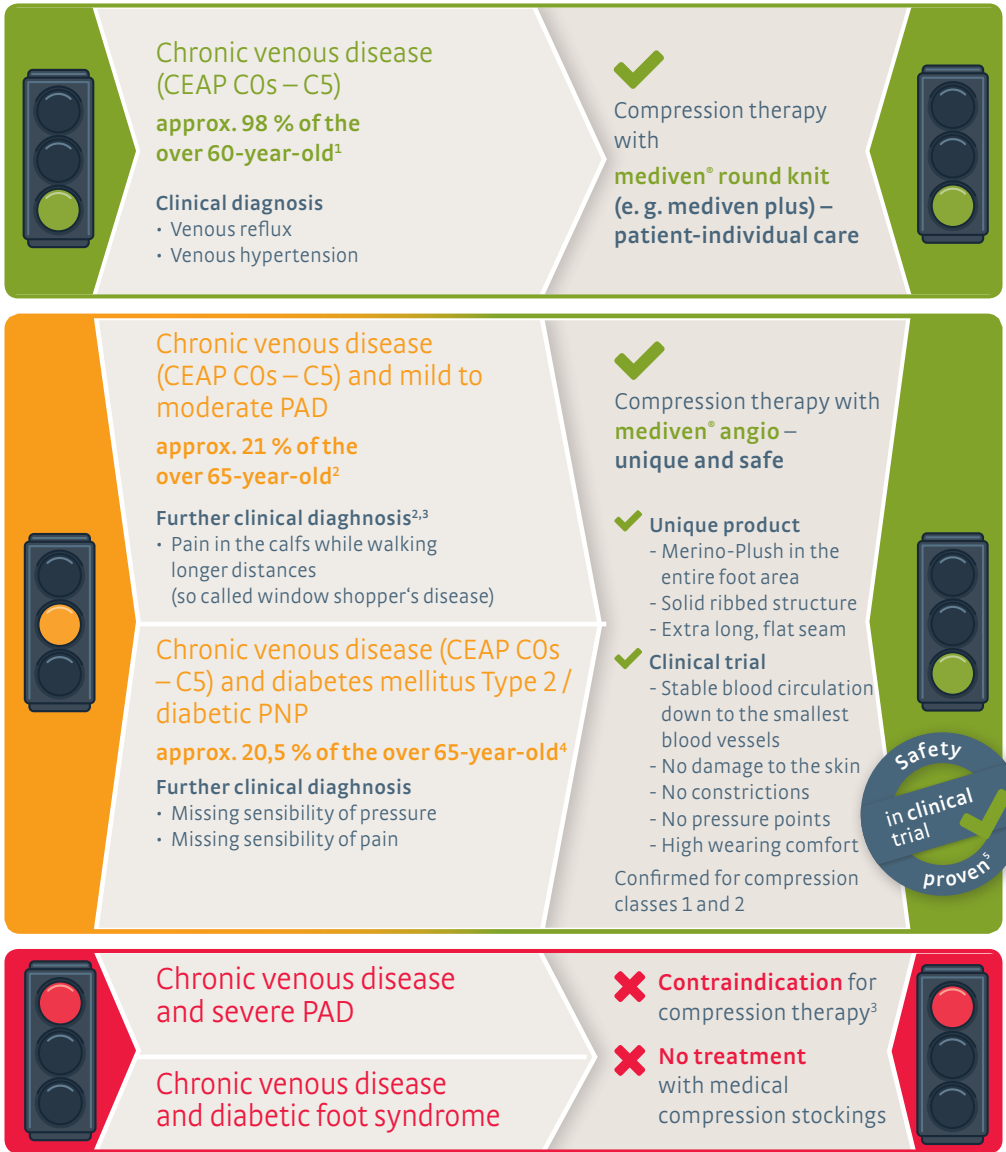

¹ Rabe E et al. Bonner Venenstudie der Deutschen Gesellschaft für Phlebologie. Phlebologie 2003;32:1–14.

² Lawall H et al. S3-Leitlinie zur Diagnostik, Therapie und Nachsorge der peripheren arteriellen Verschlusskrankheit. Online veröffentlicht unter: <https://www.awmf.org/leitlinien/detail/II/065-003.html> (Letzter Zugriff 09.12.2020)

³ Rabe E et al. S2k-Leitlinie: Medizinische Kompressionstherapie der Extremitäten mit Medizinischem Kompressionsstrumpf (MKS), Phlebologischem Kompressionsverband (PKV) und Medizinischen adaptiven Kompressionssystemen (MAK). Online veröffentlicht unter: <https://www.awmf.org/leitlinien/detail/II/037-005.html> (Letzter Zugriff 09.12.2020).





⁴ Bericht der nationalen Diabetes-Surveillance 2019. Diabetes in Deutschland. Robert Koch-Institut.

⁵ Rother U et al. Safety of medical compression stockings in patients with diabetes mellitus or peripheral arterial disease. BMJ Open Diab Res Care 2020;8:e001316.


Step 2: Consider patient-specific criteria


Recommendation:

- The higher the weight
 - The stronger the oedema (venous oedema)
 - The more severe the disease
 - The softer the connective tissue
- The stronger the stocking material should be.
 - Already meeting a single criterion requires a stronger material.
 - Please note! For lymphoedema or lipoedema, use mediven flat knit.

Material Criteria	Light (e.g., mediven elegance)	Medium (e.g., mediven plus)	Strong (e.g., mediven forte)
1. Weight			
2. Oedema (venous oedema)	—	+	++
3. Severity of the disease			






Product	mediven elegance All you see is beauty
Material	light
CCL	1, 2
Lengths	Standard: AD, AG, AT, AT / U
Standard colours	Black, Anthracite, Navy, Cashmere, Caramel, Rosé, Bronze, Beige, medi Magenta*, White*
Sizes	I – VII
Standard /made-to-measure	Standard and made-to-measure
Foot part	Open and closed toe, universal foot size, shoe sizes 36 – 43
Topband / attachment	Silicone dot topband (NoB), Silicone dot topband motif (NoM), Plain silicone topband rose (FBR), Silicone micro dot topband sensitive (MBs). All topbands are standard styles.
Special versions	Possible
Key words	Transparent, Soft Elastic foot tip, wide and soft cuff, smaller heel for wrinkle-free fit.
Care instructions	
Material composition	Elastane, Polyamide

* Delivery period up to 10 working days.


Measuring chart

mediven elegance®

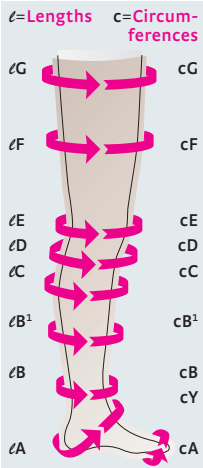Standard size seamless		I	II	III	IV	V	VI	VII	
AG with topband normal		cG	43 – 48	45 – 52	49 – 56	53 – 60	56 – 64	60 – 68	64 – 72
AG with topband wide		cG	49 – 57	53 – 62	57 – 67	61 – 72	65 – 77	69 – 81	73 – 85
 <p>ℓ=Lengths c=Circumferences</p>	cG	43 – 57	45 – 62	49 – 67	53 – 72	56 – 77	60 – 81	64 – 85	
	cF	39 – 52	41 – 56	44 – 60	47 – 65	50 – 69	53 – 73	57 – 78	
	cE	30 – 37	33 – 40	35 – 43	37 – 45	39 – 48	41 – 51	42 – 53	
	cD	27 – 33	29 – 36	32 – 39	34 – 42	36 – 45	38 – 48	40 – 50	
	cC	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51	
	cB¹	23 – 27	24 – 29	26 – 32	29 – 35	31 – 37	33 – 39	35 – 41	
	cB	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32	
	cY	26 – 31	28 – 33	29 – 35	31 – 37	32 – 38	33 – 40	34 – 42	
	cA	17 – 22	19 – 24	21 – 26	23 – 29	25 – 32	27 – 34	28 – 35	
	Lengths	ℓG normal	72 – 83 cm		Foot	Highly elastic universal foot size			
ℓG petite		62 – 71 cm							
ℓD normal		39 – 44 cm							
ℓD petite		34 – 38 cm							

mediven® for men

round-knitted

Product	mediven for men The elegant compression sock for men
Material	light
CCL	1,2
Lengths	AD
Colours	Black, Anthracite, Navy, Brown
Sizes	I – VII
Standard / made-to-measure	Standard
Foot part	Closed toe, shoe sizes 39 – 46
Key words	Anatomical toe design, right / left markings, Business Rib Design, Clima Fresh, Added Lanolin*.
Care instructions	
Material composition	Elastane, Polyamide

* Natural substance used for treating dry and chapped skin, has a broad range of applications, for example in child care and baby care. The caring protective effect remains even after washing.


Measuring chart

mediven® for men

Standard size seamless		I	II	III	IV	V	VI	VII
<p>ℓ=Lengths c=Circumferences</p>	cD	27–33	29–36	32–39	34–42	36–45	38–48	40–50
	cC	28–34	30–37	33–40	35–43	37–46	39–49	41–51
	cB¹	23–27	24–29	26–32	29–35	31–37	33–39	35–41
	cB	18–20	20–22	22–24	24–26	26–28	28–30	30–32
	cY	26–31	28–33	29–35	31–37	32–38	33–40	34–42
	cA	17–22	19–24	21–26	23–29	25–32	27–34	28–35
<p>Lengths</p>	mediven for men							
	ℓD normal	39–44 cm						
	ℓD petite	34–38 cm						

mediven® comfort

round-knitted


Product	mediven comfort All you feel is softness
Material	medium
CCL	1, 2
Lengths	Standard: AD, AG, AT
Standard colours	Black, Anthracite, Navy, Cashmere, Caramel, Rosé*, Bronze*, Beige*
Sizes	I – VII
Standard / made-to-measure	Standard and made-to-measure
Foot part	Open and closed toe, universal foot size, shoe sizes 36 – 43, double-knitted comfort sole
Topband / attachment	Silicone dot topband (NoB), Silicone dot topband motif (NoM), Plain silicone topband rose (FBR), Silicone micro dot topband sensitive (MBs). All topbands are standard styles.
Special versions	Possible
Key words	Soft material, double-knitted comfort sole, added Lanolin**.
Care instructions	
Material composition	Elastane, Polyamide

* Delivery period up to 10 working days.

** Natural substance used for treating dry and chapped skin, has a broad range of applications, for example in child care and baby care.
The caring protective effect remains even after washing.


Measuring chart

mediven® comfort

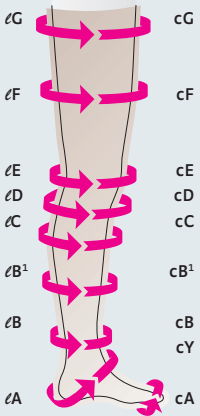Standard size seamless		I	II	III	IV	V	VI	VII	
AG with topband normal		cG	43 – 48	45 – 52	49 – 56	53 – 60	56 – 64	60 – 68	64 – 72
AG with topband wide		cG	49 – 57	53 – 62	57 – 67	61 – 72	65 – 77	69 – 81	73 – 85
<div>ℓ=Lengths c=Circumferences</div> 	cG	cG	43 – 57	45 – 62	49 – 67	53 – 72	56 – 77	60 – 81	64 – 85
	cF	cF	39 – 52	41 – 56	44 – 60	47 – 65	50 – 69	53 – 73	57 – 78
	cE	cE	30 – 37	33 – 40	35 – 43	37 – 45	39 – 48	41 – 51	42 – 53
	cD	cD	27 – 33	29 – 36	32 – 39	34 – 42	36 – 45	38 – 48	40 – 50
	cC	cC	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51
	cB1	cB1	23 – 27	24 – 29	26 – 32	29 – 35	31 – 37	33 – 39	35 – 41
	cB	cB	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32
	cY	cY	26 – 31	28 – 33	29 – 35	31 – 37	32 – 38	33 – 40	34 – 42
	cA	cA	17 – 22	19 – 24	21 – 26	23 – 29	25 – 32	27 – 34	28 – 35
	Lengths	ℓG normal	72 – 83 cm		Foot	Highly elastic universal foot size Double-knitted comfort sole			
ℓG petite		62 – 71 cm							
ℓD normal		39 – 44 cm							
ℓD petite		34 – 38 cm							

mediven plus®


round-knitted

	mediven plus All you experience is reliability
Material	medium
CCL	1, 2, 3
Lengths	Standard: AD, AG, AT, AT / U, AT / H, AT / E
Standard colours	Black, Anthracite, Navy, Cashmere, Caramel, Rosé, Bronze, Beige, medi Magenta*, White*
Sizes	I – VII
Standard / made-to-measure	Standard and made-to-measure
Foot part	Open and closed toe, universal foot size, shoe sizes 36 – 43
Topband / attachment	Silicone dot topband (NoB), Silicone dot topband motif (NoM), Plain silicone topband rose (FBR), Silicone micro dot topband sensitive (MBs). All topbands are standard styles.
Special versions	Possible
Key words	Hard-wearing and opaque, plus heel, optimised toe, ribbed sole reinforcement, 5 cm wide cuff.
Care instructions	
Material composition	Elastane, Polyamide

* Delivery period up to 10 working days.

Measuring chart

mediven plus®

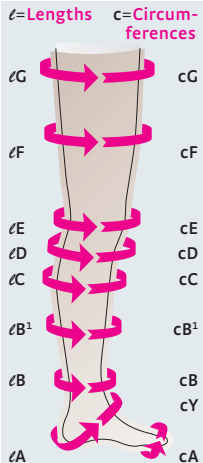Standard size seamless		I	II	III	IV	V	VI	VII	
AG with topband normal		cG	43 – 48	45 – 52	49 – 56	53 – 60	56 – 64	60 – 68	64 – 72
AG with topband wide		cG	49 – 57	53 – 62	57 – 67	61 – 72	65 – 77	69 – 81	73 – 85
<div>ℓ=Lengths c=Circumferences</div> 	cG	43 – 57	45 – 62	49 – 67	53 – 72	56 – 77	60 – 81	64 – 85	
	cF	39 – 52	41 – 56	44 – 60	47 – 65	50 – 69	53 – 73	57 – 78	
	cE	30 – 37	33 – 40	35 – 43	37 – 45	39 – 48	41 – 51	42 – 53	
	cD	27 – 33	29 – 36	32 – 39	34 – 42	36 – 45	38 – 48	40 – 50	
	cC	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51	
	cB¹	23 – 27	24 – 29	26 – 32	29 – 35	31 – 37	33 – 39	35 – 41	
	cB	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32	
	cY	26 – 31	28 – 33	29 – 35	31 – 37	32 – 38	33 – 40	34 – 42	
	cA	17 – 22	19 – 24	21 – 26	23 – 29	25 – 32	27 – 34	28 – 35	
Lengths		ℓG normal	72 – 83 cm						
		ℓG petite	62 – 71 cm						
		ℓD normal	39 – 44 cm						
		ℓD petite	34 – 38 cm						
Foot		Highly elastic universal foot size for mediven plus							

Maxi-panty top


AT mediven plus

We recommend a maxi-panty top above an abdominal circumference of

Size I – III	T = from 95 cm	H = from 110 cm
Size IV – V	T = from 120 cm	H = from 135 cm
Size VI – VII	T = from 140 cm	H = from 145 cm

mediven active®

round-knitted

Product	mediven active The hard wearing compression sock for men
Material	strong
CCL	1,2
Lengths	AD
Colours	Black, Anthracite, Navy, Brown, White
Sizes	I – VII
Standard / made-to-measure	Standard
Foot part	Closed toe, shoe sizes 38 – 48, comfort sole, seam you can't feel at the toe
Key words	Fine-ribbed design, comfort sole with reinforcement, soft and wide cuff, very robust knitted fabric.
Care instructions	
Material composition	Elastane, Polyamide

Measuring chart


mediven active®

Standard size seamless		I	II	III	IV	V	VI	VII
<p>ℓ=Lengths c=Circumferences</p>	cD	27–33	29–36	32–39	34–42	36–45	38–48	40–50
	cC	28–34	30–37	33–40	35–43	37–46	39–49	41–51
	cB¹	23–27	24–29	26–32	29–35	31–37	33–39	35–41
	cB	18–20	20–22	22–24	24–26	26–28	28–30	30–32
	cY	26–31	28–33	29–35	31–37	32–38	33–40	34–42
	cA	17–22	19–24	21–26	23–29	25–32	27–34	28–35

Lengths	mediven active	
	ℓD normal	43–49 cm
	ℓD petite	37–42 cm

mediven® forte


round-knitted

	mediven forte All you feel is protection	
Material	strong	
CCL	2, 3	
Lengths	AD, AG, AT	
Standard colours	Black, Anthracite, Navy, Cashmere, Caramel, Rosé*, Bronze*, Beige*	
Sizes	I – VII	
Standard / made-to-measure	Standard and made-to-measure	
Foot part	Open and closed toe	Small foot up to 23 cm $\hat{=}$ shoe sizes 36 – 38 Large foot from 23 cm $\hat{=}$ shoe sizes 38 – 43
Topband / attachment	Silicone dot topband (NoB), Silicone dot topband motif (NoM), Plain silicone topband rose (FBR), Silicone micro dot topband sensitive (MBs). All topbands are standard styles.	
Special versions	Possible	
Key words	Extremely hard-wearing, seamless material, opaque, Clima Fresh over the whole leg, stable soft cuff at the AD border.	
Care instructions		
Material composition	Elastane, Polyamide	

*Delivery period up to 10 working days.



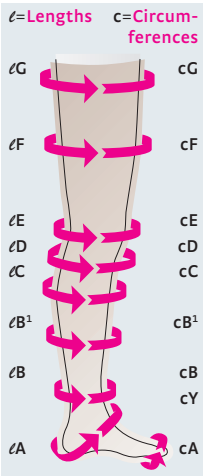


Measuring chart

mediven® forte

Standard size seamless		I	II	III	IV	V	VI	VII	
AG with topband normal		cG	44 – 47	46 – 50	50 – 54	53 – 58	56 – 61	60 – 65	64 – 70
AG with topband wide		cG	48 – 53	51 – 57	55 – 61	59 – 65	62 – 70	66 – 74	71 – 80
<div>ℓ=Lengths c=Circumferences</div> 	cG	44 – 53	46 – 57	50 – 61	53 – 65	56 – 70	60 – 74	64 – 80	
	cF	40 – 47	42 – 51	45 – 55	48 – 58	51 – 62	54 – 67	58 – 72	
	cE	30 – 37	33 – 40	36 – 42	38 – 44	40 – 46	42 – 48	44 – 50	
	cD	28 – 33	30 – 35	32 – 37	34 – 39	36 – 41	38 – 43	40 – 45	
	cC	29 – 34	31 – 36	33 – 38	35 – 40	37 – 42	39 – 44	41 – 46	
	cB¹	24 – 26	26 – 28	28 – 31	30 – 33	32 – 35	34 – 37	36 – 39	
	cB	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32	
	cY	28 – 30	29 – 32	31 – 34	33 – 36	34 – 38	36 – 40	37 – 41	
	cA	18 – 21	20 – 23	22 – 25	24 – 27	26 – 29	28 – 31	30 – 33	
	<div>Lengths</div>	ℓG normal	72 – 83 cm		<div>Foot</div> <div>Small foot up to 23 cm ≙ shoe sizes 36 – 38</div> <div>Large foot from 23 cm ≙ shoe sizes 38 – 43</div>				
ℓG petite		62 – 71 cm							
ℓD normal		38 – 43 cm							
ℓD petite		33 – 37 cm							

mediven® angio


round-knitted

Produkt	mediven angio Safety in treatment – down to the toes.
Material	medium
CCL	1, 2
Lengths (Standard & MTM*)	AD
Standard colours	Black, Caramel
Sizes	I – VII
Foot part	Closed toe: small foot up to 25 cm $\hat{=}$ shoe size 36 – 40 large foot from 25 cm $\hat{=}$ shoe size 40 – 46
Topband	Tricot border
Special versions	Not possible
Key words	Suitable for concomitant mild to moderate peripheral artery disease (PAD) and / or diabetes mellitus, thanks to Merino-Plush-Padding in the entire foot area up to above the ankles, extra long and flat toe seam, solid ribbed structure, 90°-heel, extra wide cuff, safety proven in clinical trial
Care instructions	    
Material composition	Elastane, Polyamide, Merino wool

* Delivery period up to 10 working days.

Measuring chart

mediven® angio

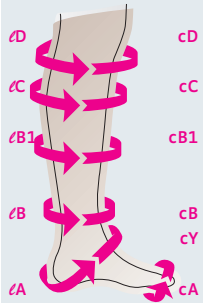Standard sizes		I	II	III	IV	V	VI	VII
<p>ℓ = lengths c = circumferences</p> 	cD	27 – 33	29 – 36	32 – 39	34 – 42	36 – 45	38 – 48	40 – 50
	cC	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51
	cB1	23 – 27	24 – 29	26 – 32	29 – 35	31 – 37	33 – 39	35 – 41
	cB	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32
	cY	26 – 31	28 – 33	29 – 35	31 – 37	32 – 38	33 – 40	34 – 42
	cA	17 – 22	19 – 24	21 – 26	23 – 29	25 – 32	27 – 34	28 – 35

Lengths	Product lengths	
	ℓD regular	39 – 44 cm
	ℓD petite	34 – 38 cm

Foot lengths	
small foot	closed toe, up to 25 cm ≙ shoe sizes 36 – 40
large foot	closed toe, from 25 cm ≙ shoe sizes 40 – 46

mediven® gusset variants

standard

Standard Gusset	mediven elegance	mediven comfort	mediven plus	mediven forte
Comfort gusset	Sizes I – VII cT 60 – 125 cH 80 – 145	Sizes I – VII cT 60 – 125 cH 80 – 135		
Crotch gusset Small abdomen and small bottom For women with a less pronounced abdomen and bottom			Sizes I – III cT 60 – 95 cH 80 – 110	Sizes I – III cT 60 – 95 cH 80 – 110
Narrow gusset at the back Small abdomen and large bottom For women with a pronounced bottom, but a relatively flat abdomen			Sizes IV – VII cT 80 – 125 cH 100 – 145 Maxi panty top Sizes I – III cT 95 – 140 cH 110 – 150	Sizes IV – VII cT 80 – 125 cH 100 – 145
Wide gusset at the back Small abdomen and very large bottom For women with a pronounced bottom but a comparatively flat abdomen			Maxi panty top Sizes IV – V cT 120 – 150 cH 135 – 160	
Extra wide gusset at the back Very large abdomen and very large bottom For women with a pronounced bottom and abdomen			Maxi panty top Sizes VI – VII cT 140 – 165 cH 145 – 170	
Maternity panty AT / U, closed Growing abdomen Suitable for pregnant women with a large, growing abdominal circumference	Sizes I – VII cT from 60 cH from 80		Sizes I – VII cT from 60 cH from 80	
Men's leotard, horizontal fly For men (with fly)			Sizes I – VII cT 60 – 125 cH 80 – 145	



mediven® gusset variants

Unless otherwise specified, the special is available in the corresponding standard sizes and made-to-measure

Gusset names	Figure descriptions
Narrow gusset back	Pronounced bottom and relatively flat abdomen <div> <div>mediven plus Size I – III</div> <div>mediven forte Size I – III</div> </div>
Narrow gusset front	Pronounced abdomen and relatively flat bottom
Wide gusset back	Very pronounced bottom and relatively flat abdomen
Wide gusset front	Very pronounced abdomen and relatively flat bottom
Variable gusset back	Narrow waist, very pronounced bottom and relatively flat abdomen
Variable gusset front	Narrow waist, very pronounced abdomen and relatively flat bottom
Continuous gusset narrow	Tenderness to pressure and slight bottom and abdomen
Continuous gusset front narrow and back wide	Large abdominal circumference and particularly pronounced bottom
Continuous gusset front wide and back narrow	Overall large abdominal circumference and a pronounced bottom
Continuous gusset wide	Large abdominal and hip circumferences
Maternity pantyhose AT / U, open	Suitable for pregnant women with large and growing abdominal circumferences and tenderness to pressure
Maternity pantyhose AT / U, closed	Suitable for pregnant women with large and growing abdominal circumferences
Men's-gusset with horizontal fly	Suitable for men (with fly)


mediven® ulcer kit

Rundstrick

Product	mediven ulcer kit Ulcer compression stocking system with 2 x 20 mmHg.
Material	Two components
CCL	2 x 20 mmHg
Lengths	Standard: AD
Colours	Liner stocking: white Outer stocking: caramel
Sizes	I – VII
Foot part	Liner stocking: closed toe , shoe sizes 36 –43 Outer stocking: closed or open toe , shoe sizes 36 – 43
Topband	Liner stocking: flat topband Outer stocking : wide topband
Special versions	Not possible
Key words	Soft, antibacterial liner stocking, easy donning sewn-in instep area, colour marking at the heel
Care instructions	Liner stocking:  Outer stocking: 
Material composition	Elastane, Polyamide

Measuring chart

mediven® ulcer kit

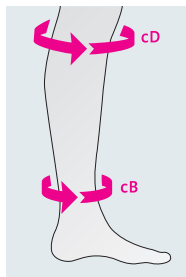

Measurements / Size table		Size I	Size II	Size III	Size IV	Size V	Size VI	Size VII
Leg measurements (in cm)								
Ankle cB		18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32
Below knee cD		27 – 33	29 – 36	32 – 39	34 – 42	36 – 45	38 – 48	40 – 50
Colour marking	Heel	pink	blue	yellow	red	grey	green	brown
Length (in cm)		petite		regular				
Ø		34 – 38		39 – 44				

Your order includes

- 2 mediven ulcer liner stockings
- 1 mediven ulcer plus outer stocking
- 1 donning aid for stockings with open toe (only when ordering mediven ulcer kit open toe)



mediven harmony®

round-knitted

Product	mediven harmony Round-knitted armsleeve, flat-knitted hand part
CCL	1, 2
Lengths	Standard length: AG, CG (regular and extra wide) Short length: AG, CG with topband (regular and extra wide) AC ¹ hand piece or glove
Colours	Black*, Caramel, Sand
Sizes	I – VIII, armsleeve II – VII, hand part / fingers
Standard / made-to-measure	Standard
Topband / attachment	Narrow napped silicone topband, shoulder cap
Key words	Comfortable and effective, double ellipses on the point and in the crease of the elbow
Care instructions	    
Material composition	Elastane, Polyamide

*Black: minimum order quantity of 20 pieces per size and variant.

Armsleeve	Size I		Size II		Size III	
Circumference (cm)	Regular	Extra wide	Regular	Extra wide	Regular	Extra wide
cG	25,5 – 28,0	28,5 – 31,0	28,0 – 30,5	30,5 – 33,5	30,0 – 33,0	33,0 – 36,0
cF	23,5 – 25,5	25,5 – 28,0	25,0 – 27,5	27,5 – 30,0	27,5 – 30,0	30,0 – 33,0
cE	23,0 – 25,0	24,0 – 26,0	24,5 – 26,5	26,5 – 29,0	26,0 – 28,5	28,5 – 31,0
cD	22,0 – 24,0	24,0 – 26,0	23,0 – 25,0	25,0 – 27,5	24,0 – 26,0	26,0 – 28,5
cC 1 = 100% **	15,0 – 16,5	15,0 – 16,5	16,5 – 18,0	16,5 – 18,0	18,0 – 19,5	18,0 – 19,5
cC	13,5 – 14,5	13,5 – 14,5	14,5 – 15,5	14,5 – 15,5	16,0 – 17,5	16,0 – 17,5
cA	15,0 – 18,0	15,0 – 18,0	16,0 – 19,0	16,0 – 19,0	17,0 – 20,0	17,0 – 20,0
Armsleeve	Size IV					
Circumference (cm)	Regular	Extra wide				
cG	32,0 – 35,0	35,0 – 38,5				
cF	29,5 – 32,0	32,0 – 35,0				
cE	27,5 – 30,0	30,0 – 33,0				
cD	25,0 – 27,5	28,0 – 30,5				
cC 1 = 100% **	19,5 – 21,0	19,5 – 21,0				
cC	17,5 – 19,0	17,5 – 19,0				
cA	19,0 – 22,0	19,0 – 22,0				

Armsleeve	Size V		Size VI		Size VII	
Circumference (cm)	Regular	Extra wide	Regular	Extra wide	Regular	Extra wide
cG	33,5 – 36,5	37,0 – 40,5	36,0 – 39,5	39,5 – 43,0	38,0 – 41,5	41,5 – 45,5
cF	31,0 – 34,0	34,0 – 37,0	33,0 – 36,0	36,0 – 39,5	35,0 – 38,5	38,5 – 42,0
cE	29,0 – 31,5	31,5 – 34,5	30,5 – 33,5	33,5 – 36,5	32,5 – 35,5	35,5 – 39,0
cD	27,0 – 29,5	29,5 – 32,0	28,5 – 31,0	31,0 – 34,0	30,5 – 33,5	33,5 – 36,5
cC 1 = 100% **	21,0 – 23,0	21,0 – 23,0	23,0 – 25,0	23,0 – 25,0	25,0 – 27,5	25,0 – 27,5
cC	19,0 – 20,5	19,0 – 20,5	20,5 – 22,5	20,5 – 22,5	22,5 – 24,5	22,5 – 24,5
cA	20,0 – 23,0	20,0 – 23,0	22,0 – 25,0	22,0 – 25,0	24,0 – 27,0	24,0 – 27,0

Armsleeve	Size VIII***	
Circumference (cm)	Regular	Extra wide
cG	40,5 – 44,5	45,0 – 49,5
cF	37,5 – 41,0	41,5 – 45,5
cE	35,0 – 38,5	38,5 – 42,0
cD	33,0 – 36,0	36,0 – 39,5
cC 1 = 100% **	27,5 – 30,0	27,5 – 30,0
cC	25,0 – 27,5	25,0 – 27,5
cA	26,0 – 29,0	26,0 – 29,0



Length (cm)	Standard length	Short length
	41 – 47	36 – 40

** Measurement CC1 is 5 cm above CC

*** Size VIII: Minimum order quantity of 20 pieces per colour and variant;
delivery period approx. 4 weeks

Measuring chart

hand part

Variants	
Seamless AC1 hand piece without fingers, with thumb base	Seamless AC1 glove with compressive fingers
	

Length / Hand part (cm), all sizes

AB	3,5 – 4,5
AC	10,0 – 12,0
AC1	15,0 – 16,0

Circumference / Hand part (cm)						
Sizes	II	III	IV	V	VI	VII
cC 1 = 100 % *	15,0 – 18,0	17,0 – 20,0	19,0 – 21,5	21,0 – 23,5	23,0 – 26,0	25,5 – 30,0
cC	13,5 – 15,5	15,0 – 18,0	17,0 – 20,0	19,0 – 21,0	20,0 – 23,5	22,5 – 27,5
cB	15,5 – 20,0	17,5 – 21,0	19,5 – 23,0	20,5 – 24,0	22,5 – 26,0	24,5 – 30,0
cA	15,0 – 19,0	17,0 – 20,0	19,0 – 22,0	20,0 – 23,0	22,0 – 25,0	24,0 – 29,0



Length / Fingers (cm)

Thumb	3,5 – 4,5
Index finger	4,5 – 5,5
Middle finger	5,5 – 6,5
Ring finger	4,5 – 5,5
Little finger	3,5 – 4,5


* Measurement CC1 is 5 cm above CC

Circumference / Fingers (cm) / Sizes

Sizes		II	III	IV	V	VI	VII
Thumb	Bottom cX	7,0	7,5	8,0	8,5	9,0	10,0
	Top cZ	6,0	6,5	6,5	7,0	7,5	8,0
Index finger	Bottom cX	6,5	7,0	7,5	8,0	8,5	9,5
	Top cZ	5,0	5,5	6,0	6,0	6,5	7,0
Middle finger	Bottom cX	6,5	7,0	7,5	8,0	8,5	9,5
	Top cZ	5,0	5,5	6,0	6,0	6,5	7,0
Ring finger	Bottom cX	6,0	6,5	7,0	7,0	7,5	9,0
	Top cZ	5,0	5,0	5,5	5,5	6,0	6,5
Little finger	Bottom cX	5,5	6,0	6,0	6,5	7,0	8,0
	Top cZ	4,5	4,5	5,0	5,0	5,5	6,0

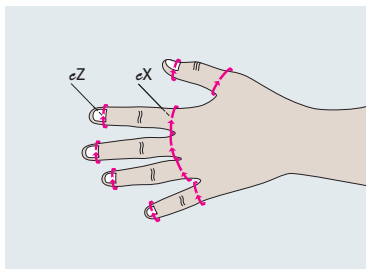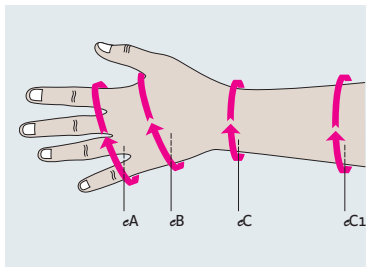

Circumference eZ is measured at the nailbed.
Circumference eX is measured at the base of the finger.

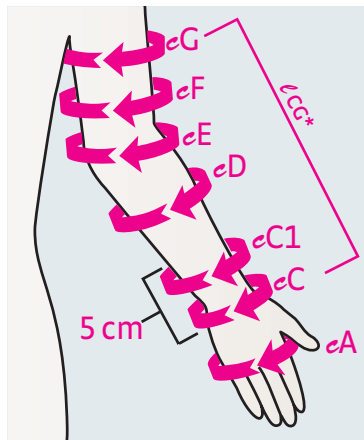

Measurement median harmony

* Length eCG is measured on the inside of the arm from the wrist to the anterior axillary fold.

mediven® flat knit care compass leg


medi tip:

We refer to the current German S2k guideline Medical Compression Therapy¹, the S1 guideline Lipoedema² and the S2k guideline Diagnostics and Therapy of Lymphoedema³. For the supply of compression therapy aids, the relevant valid guidelines must be observed.


Step 1: Indication

Indication	mediven® 550 leg	mediven® cosy	mediven mondi®	circaid® juxtafit®
Primary and secondary lymphoedema	●	●	●	●
Lipoedema	●	●	●	
Scar treatment	●	●	●	
Post-operative and post-traumatic oedema	●	●	●	
Mixed forms such as phlebo-lymphoedema, lipo-lymphoedema	●	●		●

Step 2: Consider the stage


















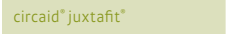


Indication	Diagnosis	Characteristics	mediven® 550 leg	mediven® cosy	mediven mondi®	circaid® juxtafit®
Lipoedema 	Stage 1	Smooth skin surface, uniform thickening of the subcutis at certain points	●	●	●	
	Stage 2	Wave-like skin surface, knot-like structures in thickened subcutis	●	●	●	
	Stage 3	Coarse nodular skin surface, circumferential augmentation with overhanging tissue parts	●	●		

Step 2: Consider the stage

Indication	Diagnosis	Characteristics	mediven® 550 leg	mediven® cosy	mediven mondi®	circaid® juxtafit®
Lymphoedema 	Stage 1	Spontaneously reversible, oedema of soft consistency; elevation reduces swelling, pitting of the skin	●	●	●	●
	Stage 2	Spontaneously irreversible, hardened oedema; elevation does not eliminate swelling, skin does not pit	●	●	●	●
	Stage 3	Unformed swelling of extreme dimensions, severely hardened skin, formation of blisters and fistulae	●			●

Because of the material, the fit and the product features, we recommend the mediven flat knit products for patients with the following characteristics*.

Step 3: Consider patient-specific criteria

Criteria	Material*	Light e.g. mediven mondi®	Medium e.g. mediven® cosy	Strong e.g. mediven® 550 leg
				
Connective tissue		normal to firm	soft to firm	very soft to hardened
Weight				
Deepened skin wrinkles / calibre jumps		none	moderate	strong
Severity of the disease				
Recommended mediven® compression stockings:				
				
				
Recommended for patients who are not able to apply compression stockings:				

* wall stability

Round-knitted

Flat-knitted

Recommendation:

- The **softer** the connective tissue,
- The **higher** the weight,
- The **more** deepened skin wrinkles,
- The **more** severe the disease.

→ The **stronger** the stocking material should be.

→ The most prominent condition indicates stocking choice.

Leg garment	Oedema stage	Knitting technique	S = Standard M = Made-to-measure
mediven 550 leg	1, 2, 3	flat knit	S, M
mediven cosy	1, 2, 3	flat knit	M
mediven mondi	1, 2	flat knit	S, M
circaid juxtafit premium and essentials leg	1, 2, 3	adjustable multi components textile	S, M

¹ Rabe E et al. S-k guideline: Medical compression therapy of the extremities with medical compression stockings (MCS), phlebological compression bandages (PKV) and Medical Adaptive Compression Systems (MAK). Published online at: <https://www.awmf.org/leitlinien/detail/II/037-005.tml> (Last accessed 18.03.2020)

² Reich-Schupke S et al. S1 guideline lipoedema. Published online at: <https://www.awmf.org/leitlinien/detail/II/037-012.html> (Last accessed 18.03.2020)

³ Oberlin M et al. S2k guideline: Diagnosis and therapy of lymphoedema. Published online at: <https://www.awmf.org/leitlinien/detail/II/058-001.html> (Last accessed 18.03.2020)

⁴ List of aids of the GKV-Spitzenverband, PG 17, Link: https://hilfsmittel.gkv-spitzenverband.de/produktartlisteZurPG_input.action?paramGruppel=17 (Last accessed 18.03.2020)

*Please note that this recommendation does not replace the diagnosis of a doctor or the advice of a medical specialist. The choice of the optimal compression stocking depends on the indication and the individual needs of the patient.

mediven® flat knit care compass arm

medi tip:


We refer to the current German S2k guideline Medical Compression Therapy Compression Therapy¹, the S1 guideline Lipoedema² and the S2k guideline Diagnostics and Therapy Lymphoedema³.

For the supply of compression therapy aids, the relevant valid guidelines must be observed.


Step 1: Indication

Indication	mediven® 550 arm	mediven® mondi esprit	circaid® juxtafit®
Primary and secondary Lymphoedema	●	●	●
Lipoedema	●	●	
Scar treatment	●	●	
Post-operative and post-traumatic oedema	●	●	

Step 2: Consider stage

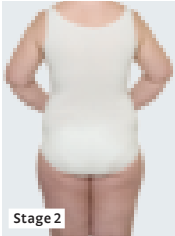








Indication	Diagnosis	Characteristics	mediven® 550 arm	mediven® mondi esprit	circaid® juxtafit®
	Stage 1	Smooth skin surface, uniform thickening of the subcutis at certain points	●	●	
	Stage 2	Wave-like skin surface, knot-like structures in thickened subcutis	●	●	
	Stage 3	Coarse nodular skin surface, circumferential augmentation with overhanging tissue parts	●		

Step 2: Consider stage

Indication	Diagnosis	Characteristics	mediven® 550 arm	mediven® mondi esprit	circaid® juxtafit®
 <p>Stage 3</p>	Stage 1	Spontaneously reversible, oedema of soft consistency; elevation reduces swelling, pitting of the skin	●	●	●
	Stage 2	Spontaneously irreversible, hardened oedema; elevation does not eliminate swelling, skin does not pit	●	●	●
	Stage 3	Unformed swelling of extreme dimensions, severely hardened skin, formation of blisters and fistulae	●		●

Because of the material, the fit and the product features, we recommend the mediven flat knit products for patients with the following characteristics*.

Step 3: Consider patient-specific criteria

Criteria	Material*	Light / Medium e.g. mediven® mondi esprit	Strong e.g. mediven® 550 arm
			
Connective tissue		normal to firm	very soft to hardened
Weight			
Deepened skin wrinkles / calibre jumps		none / moderate	moderate / strong
Severity of the disease			
Recommended mediven® compression armsleeves:		mediven® mondi esprit	
			mediven® 550 arm
Recommended for patients who are not able to apply compression armsleeves:		circaid® juxtafit®	

* wall stability

Recommendation:

- The **softer** the connective tissue,
- The **higher** the weight,
- The **more** deepened skin wrinkles,
- The **more** severe the disease.

→ The **stronger** the stocking material should be.

→ The most prominent condition indicates stocking choice.

Arm garment	Oedema stage	Knitting technique	S = Standard M = Made-to-measure
mediven 550 arm	1, 2, 3	flat knit	M
mediven esprit	1, 2	flat knit	S, M
mediven harmony	1	round knit	S
circaid juxtafit essentials arm	1, 2, 3	adjustable multi-component textile	S, M

¹ Rabe E et al. S+k-Leitlinie: Medizinische Kompressionstherapie der Extremitäten mit Medizinischem Kompressionsstrumpf (MKS), Phlebologischem Kompressionsverband (PKV) und Medizinischen adaptiven Kompressionssystemen (MAK). Online veröffentlicht unter: <https://www.awmf.org/leitlinien/detail/II/037-005.tml> (Letzter Zugriff 18.03.2020)

² Reich-Schupke S et al. S1-Leitlinie Lipödem. Online veröffentlicht unter: <https://www.awmf.org/leitlinien/detail/II/037-012.html> (Letzter Zugriff 18.03.2020)






³ Oberlin M et al. S2k-Leitlinie: Diagnostik und Therapie der Lymphödeme. Online veröffentlicht unter: <https://www.awmf.org/leitlinien/detail/II/058-001.html> (Letzter Zugriff 18.03.2020)

⁴ Hilfsmittelverzeichnis des GKV-Spitzenverband, PG 17, Link: https://hilfsmittel.gkv-spitzenverband.de/produktartlisteZurPG_input.action?paramGruppeld=17 (Letzter Zugriff 18.03.2020)

* Please note that this recommendation does not replace the diagnosis of a doctor or the advice of a medical specialist. The choice of the optimal compression stocking depends on the indication and the individual needs of the patient.

mediven® 550 leg


flat-knitted

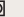Product	mediven 550 leg* Strong. Very high wall stability. Reliable hold.
CCL	2
Lengths	2 lengths for AD / calf and AD / calf petite
Standard colour	Caramel
Sizes	IV – XIV
Standard / made-to-measure	Standard and made-to-measure
Foot part	Open toe straight ending
Topband /attachment	Silicone dot topband
Key words	Clima Comfort, Clima Fresh
Care instructions	    
Material composition	Elastane, Polyamide

* Made-to-measure products and variants please see our product catalogue.

mediven mondi®

flat-knitted

Product	mediven mondi* Soft. Fine. Looks great.
CCL	2
Lengths	2 lengths for AD / calf and AG / thigh
Standard colour	Caramel
Sizes	I – VII
Standard /made-to-measure	Standard and made-to-measure
Foot part	Open toe straight ending
Topband /attachment	Silicone dot topband
Key words	Clima Comfort, Clima Fresh
Care instructions	
Material composition	Elastane, Polyamide

* Made-to-measure products and variants please see our product catalogue.

Measuring chart

mediven mondi®

Standard size		I	II	III	IV	V	VI	VII
<p>ℓ=Lengths c=Circumferences</p>	cG with topband	44 – 47	46 – 50	50 – 54	53 – 58	56 – 61	60 – 65	64 – 70
	cG	44 – 53	46 – 57	50 – 61	53 – 65	56 – 70	60 – 74	64 – 80
	cF	40 – 47	42 – 51	45 – 55	48 – 58	51 – 62	54 – 67	58 – 72
	cE	30 – 37	33 – 40	36 – 42	38 – 44	40 – 46	42 – 48	44 – 50
	cD	28 – 33	30 – 35	32 – 37	34 – 39	36 – 41	38 – 43	40 – 45
	cC	29 – 34	31 – 36	33 – 38	35 – 40	37 – 42	39 – 44	41 – 46
	cB ₁	24 – 26	26 – 28	28 – 31	30 – 33	32 – 35	34 – 37	36 – 39
	cB	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32
	cY	28 – 30	29 – 32	31 – 34	33 – 36	34 – 38	36 – 40	37 – 41
	cA	18 – 21	20 – 23	22 – 25	24 – 27	26 – 29	28 – 31	30 – 33
	ℓG standard	72 – 83 cm						
	ℓG short	62 – 71 cm						
	ℓD standard	38 – 43 cm						
	ℓD short	33 – 37 cm						

mediven esprit®

flat-knitted

Product	mediven esprit* Soft. Fine. Looks great.
CCL	1, 2
Lengths	Standard length: AG, CG Long length: AG, CG AC ¹ hand piece
Standard colour	Caramel
Sizes	I – V
Standard / made-to-measure	Standard and made-to-measure
Topband / attachment	Narrow silicone dot topband
Key words	Clima Comfort, Clima Fresh
Care instructions	
Material composition	Elastane, Polyamide

* Made-to-measure products and variants please see our product catalogue.

Measuring chart

mediven esprit®

Variants

Styles


CG / SN

CG / SN long

AG / SN

AG / SN long

AC₁

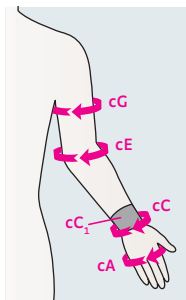

c = circumferences

Standard sizes

Circumferences in cm	Sizes	I	II	III	IV	V
	cA	17.0 – 19.0	18.0 – 20.0	19.0 – 22.0	21.0 – 24.0	23.0 – 27.0
	cC	15.0 – 17.0	17.0 – 19.0	19.0 – 21.0	21.0 – 23.0	23.0 – 25.0
	cE	25.0 – 27.5	27.0 – 29.5	29.0 – 32.0	31.0 – 34.0	33.0 – 37.0
	cG	28.0 – 32.0	30.0 – 34.0	32.0 – 37.0	35.0 – 40.0	38.0 – 44.0

Lengths

CG standard 40 – 44 cm

CG long 44 – 48 cm

Round-knitted

Flat-knitted

Measuring chart

mediven esprit®

Circumference / Hand part (cm)					
Sizes	I	II	III	IV	V
a	17 – 20	19 – 22	21 – 24	23 – 26	25 – 28
b	18 – 21	20 – 23	22 – 25	24 – 27	23 – 25
c	15 – 17	17 – 19	19 – 21	21 – 23	23 – 25
c ₁	15.5 – 18	17.5 – 20	19.5 – 22	21.5 – 24	23 – 26

Length / Hand part (cm), all sizes	
AB	3.5 – 4.5
AC	10.0 – 12.0
AC ₁	15.0 – 16.0

Measuring chart

mediven esprit®

Circumference / Fingers (cm)						
Variant AC1		I	II	III	IV	V
Thumb	top	4.5	5.5	6.5	7.5	8.5
	bottom	6.5	7.5	8.5	9.5	10.5
Index finger	top	4.5	5.5	6.5	7.5	8.5
	bottom	5.5	6.5	7.5	8.5	9.5
Middle finger	top	4.5	5.5	6.5	7.5	8.5
	bottom	5.5	6.5	7.5	8.5	9.5
Ring finger	top	4.0	5.0	6.0	7.0	8.0
	bottom	5.0	6.0	7.0	8.0	9.0
Little finger	top	3.5	4.5	5.5	6.5	7.5
	bottom	4.5	5.5	6.5	7.5	8.5

Length / Fingers (cm)	
Thumb	3.5 – 4.5
Index finger	4.5 – 5.5
Middle finger	5.5 – 6.5
Ring finger	4.5 – 5.5
Little finger	3.5 – 4.5

mediven® standard colours

standard

Round-knitted Models	mediven elegance	mediven for men	mediven comfort	mediven plus	mediven active	mediven forte	mediven angio
Black	●	●	●	●	●	●	●
Anthracite	●	●	●	●	●	●	
Navy	●	●	●	●	●	●	
Cashmere	●		●	●		●	
Caramel	●		●	●		●	●
Rosé	●		○	●		○	
Bronze	●		○	●		○	
Beige	●		○	●		○	
medi Magenta	○			○			
White	○			○	●		
Brown		●			●		
Flat-knitted	mediven 550 leg	mediven mondi	mediven esprit	mediven harmony	Note: Made-to-measure products and variants are also available in trend colours, with Design- or Fashion-Elements and can also be embellished with crystal. For a complete overview please see our product catalogue.		
Sand	○	○	○	●			
Caramel	●	●	●	●			
Black	○	○	○	●			
Cashmere	○	○	○				
Anthracite	○	○	○				
Navy	○	○	○				

○ Delivery period up to 10 working days.

medi Butler

Donning aids

Features

- Simplifies donning of all compression stockings and pantyhose
- Prescribable
- Stability gives safety
- Simple handling
- Different special models for every purpose


medi Big (XXL) Butler

For easier donning of large stocking sizes, recommended for sizes V, VI, VII.



medi Butler with variable handles

For patients with considerably restricted freedom of movement, the handle length can be adjusted individually. Accessories: screwdriver

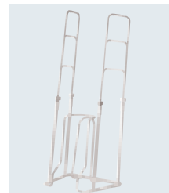


medi Arm Butler

Simplifies daily donning without assistance, frame diameter 11 cm.



medi Butler for pantyhose

For simple donning of compression pantyhose.

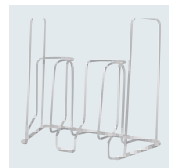

medi Donning and doffing aids


medi 2in1

Thanks to simple handling, the stocking glides up over the instep and heel almost by itself. This means far less effort is required when donning and doffing compression stockings and pantyhose.



medi slide

Easy donning of compression stockings with open toe. Thanks to the smooth material the compression stocking glides over the heel almost by itself.

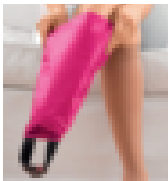

medi Butler Off

Slip out of compression stockings easily and without tiresome bending.



medi skin care

Compression stockings are the single most important form of treatment for venous disease. Skin care is a very important factor in helping to make the treatment with compression stockings as comfortable and satisfactory as possible for your patient.



medi day - cooling gel with horse chestnut.

For a fresh feeling throughout the day.

- Refreshing and cooling effect
- Well-nourished, smooth skin
- Reduces the feeling of tightness
- Supports user comfort of medical compression stockings



medi fresh – refreshing spray with menthol.

For a hint of freshness throughout the day.

- Feel refreshed straight away
- Can be used several times a day
- Gentle to the stocking fabric
- Easy to use everyday



medi night - skincare cream with regenerative witch hazel.

To soothe and rejuvenate during the night.

- Soothes the skin and also supports regeneration



medi soft – skincare with horse chestnut and moisturising urea.

Gentle foam lotion to care for the skin around the clock.

- Well nourished, smooth skin
- Reduces the feeling of tightness
- Supports user comfort of medical compression stockings
- Provides additional support for fabrics and adhesive tapes

medi textile care

Gentle care for hand washing your garments.



medi clean

For the easy care of your compression garments – anytime, anywhere.
The special washing liquid is specifically compatible for use with compression garments, cleaning with the utmost of care.



medi spot ex

The clever solution for stain-free compression garments.
Use before the wash, the spot removal stick medi spot ex removes shoe marks from the garment's toe or heel.



medi fix

Adhesive lotion for compression stockings and sleeves with tricot border for a secure fit.

Intended purposes

Intended purpose round and flat knitted medical compression stocking:

Round- or flat-knitted medical stockings are used for compression of the lower extremities, mainly for the treatment of disorders of the venous and lymphatic systems

Intended purpose medi Butler (leg):

Donning aids for medical/clinical compression socks and stockings.

Intended purpose medi Butler (arm):

Donning aid for medical/clinical compression sleeves.

Intended purpose mediven ulcer kit :

Circular-knit compression stocking system used for compression of the lower extremities, primarily in the treatment of disorders of the venous system.

mediven harmony:

Round-knitted medical compression stocking used for compression of the upper extremities, mainly for the treatment of disorders of the venous and lymphatic systems.

mediven harmony hand part:

Flat-knitted medical compression stocking used for compression of the upper extremities, mainly for the treatment of disorders of the venous and lymphatic systems.

Headquarters
medi GmbH & Co. KG
Medicusstraße 1
95448 Bayreuth
Germany
T +49 921 912-0
F +49 921 912-781
export@medi.de
www.medi.de/en



0097E17/06.2021